


RTS series

SOFTWARE

Features of RTS Series

Ties the electronic balance to windows

With this software, you can take in the RS-232C output data from SHINKO DENSHI electronic balance by the standard application software for Windows. Data from the electronic balance are transferred with the keyboard image for the applicable application.

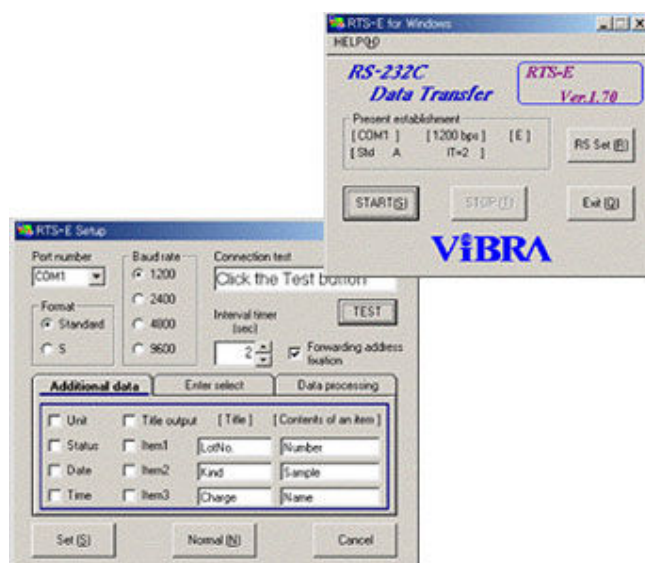

Personal Computer / Operating System

Operating system	Windows 98 Windows 98SE Windows NT4.0 (SP3 or more) Windows Me Windows 2000 Professional Windows XP Professional Edition Windows XP Home Edition
Recommended model ^(*)	installed i486DX4 66MHz CPU or more.
Memory ^(*)	16MB and more
Display of the resolution ^(*)	640×480 or more

Free space to install	more than 20MB
Other	RS-232C port (COM1-8) CD-ROM drive

- *In Windows Me/2000/XP, follow each specification recommended by OS.

RS-232C Conditions

It copes with our electronic balances.

Please set up the output control function of the electronic balance according to the condition of input data.

Standard condition

COM port	1 to 8
Baud rate	1200 / 2400 / 4800 / 9600 bps
Format	Standard / S
Parity	None / Odd / Even
Additional data	Unit / Status / Date / Time / Item 1 / Item 2 / Item 3

Fixed condition

Synchronous form	Asynchronous
Communication form	Full duplex

Protocol	Without a protocol
Start bit	1 bit
Data bit	8 bit (Standard) / 7 bit (S)
Stop bit	2 bit (Standard) / 1 bit (S)
Terminator	CR LF
